

VEVOR[®]

TOUGH TOOLS, HALF PRICE

Technical Support and E-Warranty Certificate www.vevor.com/support

RED LIGHT THERAPY BELT

MODEL: PLR-10S

We continue to be committed to provide you tools with competitive price.

"Save Half", "Half Price" or any other similar expressions used by us only represents an estimate of savings you might benefit from buying certain tools with us compared to the major top brands and does not necessarily mean to cover all categories of tools offered by us. You are kindly reminded to verify carefully when you are placing an order with us if you are actually saving half in comparison with the top major brands.

VEVOR®

TOUGH TOOLS, HALF PRICE

RED LIGHT THERAPY BELT

MODEL: PLR-10S

RED-LIGHT-THERAPY FOR HUMAN & ANIMALS!

THANK YOU FOR PURCHASING OUR LIGHT, WE HOPE YOU CAN REAP BENEFITS OF RED LIGHT THERAPY AT HOME.



NEED HELP? CONTACT US!

Have product questions? Need technical support? Please feel free to contact us:

Technical Support and E-Warranty Certificate
www.vevor.com/support

This is the original instruction, please read all manual instructions carefully before operating. VEVOR reserves a clear interpretation of our user manual. The appearance of the product shall be subject to the product you received. Please forgive us that we won't inform you again if there are any technology or software updates on our product.

Say hello to your red light!

Before starting your first session, please make sure to read this manual as it may explain important operating procedures.

Notes: Each led diode inside with three chips at 1 pcs 660nm and 2 pcs 850nm.

Near Infrared 850nm is invisible to the human naked eyes.

Specifications:

Rating power	80W
Actual Power	40W
Product Size	810*300*6 MM
LED QTY	400*0.5W
Net Weight	1.65KG
Wavelength	1*660nm;2*850nm
EMF	OUT
Irradiance	Max: 35mW/cm ²
Timer	M5/10/15/20/25/30min
Material	Fabric
Lifespan	Over 50,000 Hours

This package includes:

- | | |
|---------------------------------|---------------------------|
| 1* Red light therapy pad | 1*Mobile pocket |
| 1* Power cord | 1* Stretchy straps |
| 1* Adaptor | 1* User manual |
| 1* Controller | 1* Goggles |

The red light device uses combo red 660nm and Infrared 850nm LEDs which is chosen to give the ideal wavelengths found in the scientific research of the last three decades.

660nm of red light



This light can be seen by the human eyes and targets skin cells.

850nm of near infrared light

This light cannot be seen by the human eyes and targets deep tissue

Red Light Therapy Belt

400 pcs Red Lights + 800 pcs Infrared Lights, unique pulse function, the effect is more obvious

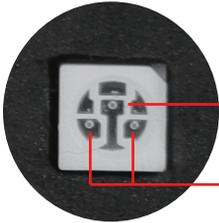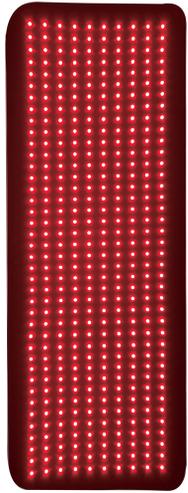

660nm Red Light

850nm Near Infrared Light

CONTROL PANEL



Red Light Energy
L0,L1,L2,L3,L4

Timer
5,10,15,20,25,30 minutes options

Red Light Energy
L4,L3,L2,L1,L0

Power /Mode Switch
Power On, Pulse Mode,Power Off

Benefits backed by science:

Research shows that the wavelengths in the spectrum of 660-900nm show optimal biotical response.

- * Muscle and Joints Pain
- * Weight loss
- * Aging spots and wrinkle removal
- * Hair Re-growth
- * Decreases pain
- * Enhanced blood circulator
- * Speed wound healing
- * Anti-inflammatory
- * Increase collagen


Even though this technology don't have side effects , we suggest taking the following precautions:

- * This device is not waterproof !
- * Not intended for the use by children and pregnant woman.
- * Do not exceed 30 Min.on one same treatment.

Application

Here show you some examples how to use the red light device.You and your pet scan enjoy the red light at different ways for like shoulder , back, feet knee ect.

LED THERAPY LIGHT PAD FOR JOINTS PAIN RELIEF

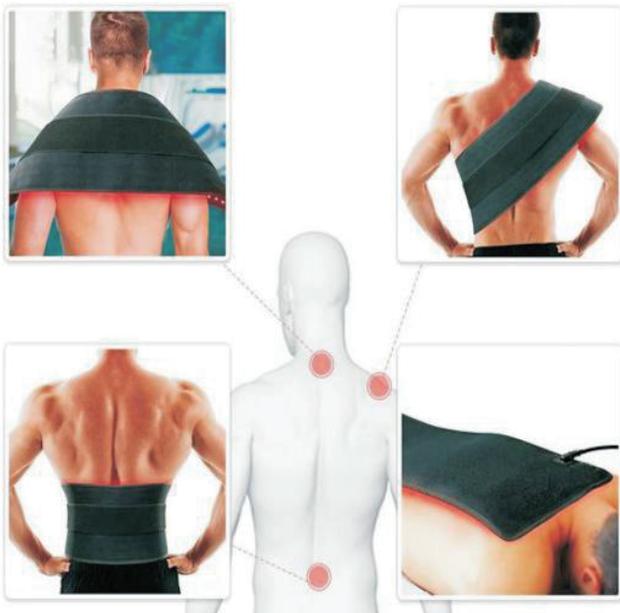



Thank you for choosing our red light therapy panel. We really appreciate your support and hope you'll enjoy this device.



SAFETY CAUTION:

1. Do not make direct view or observation to the laser beams via any optical instruments.
2. Keep out of reach of children to prevent injuries to their eyes.
3. It's forbidden for pacemaker users since the device may interfere with the operation of pacemaker.
4. During the treatment period, the diet structure should be improved and proper exercises are necessary to enhance the metabolism.
5. When using the device, please keep light emitting ports clean. Otherwise, the light losses will increase and then affect effects.
6. The symptoms like thirst, dry eyes, dizziness or rise of pulse rate after use are normal phenomena.

CORRECT DISPOSAL



This product is subject to the provision of European Directive 2012/19/EU. The symbol showing a wheelie bin crossed through indicates that the product requires separate refuse collection in the European Union. This applies to the product and all accessories marked with this symbol. Products

marked as such may not be discarded with normal domestic waste, but must be taken to a collection point for recycling electrical and electronic devices.

FCC Information:

CAUTION: Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment!

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions:

- 1) This product may cause harmful interference.
- 2) This product must accept any interference received, including interference that may cause undesired operation.

WARNING: Changes or modifications to this product not expressly approved by the party responsible for compliance could void the user's authority to operate the product.

Note: This product has been tested and found to comply with the limits for a Class B digital device pursuant to Part 15 of the FCC Rules, These limits are designed to provide reasonable protection against harmful interference in a residential installation.

This product generates, uses and can radiate radio frequency energy, and if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this product does cause harmful interference to radio or television reception, which can be determined by turning the product off and on, the user is encouraged to try to correct the interference by one or more of the following measures.

- Reorient or relocate the receiving antenna.
- Increase the distance between the product and receiver.
- Connect the product to an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for assistance.

VEVOR[®]

TOUGH TOOLS, HALF PRICE

Technical Support and E-Warranty Certificate

www.vevor.com/support